

TOSHIBA

GENESIS Edition Transforming CT

Aquilion ONE
GENESIS Edition



Transforming clinical confidence
Transforming patient experience
Transforming your workspace

GENESIS Edition Transforming CT

Brought to you by the leaders in area detector technology

Aquilion ONE™ / GENESIS Edition goes beyond the evolution of dynamic volume CT. Intensive clinically focused research and breakthrough technological developments have culminated in a CT system with industry-leading spatial resolution and reduced radiation dose requirements.

GENESIS Edition maximizes the patient experience during CT examinations, and through intelligent examination protocols, provides excellent image quality with low radiation and contrast dose tailored to each and every patient.

Intelligent technology for increased patient safety and superior patient care.





Transforming clinical confidence


The right balance between image quality and dose for every patient, from the youngest to the largest



GENESIS Edition transforms routine imaging to new levels of image detail and low contrast resolution – balanced for each clinical question at the right dose.

A completely redesigned X-ray system from photon generation to beam distribution and detection is the basis of PUREVISION Optics. This results in a better balance between image quality and dose.

Adaptive scatter correction removes scatter through intelligent modeling that preserves more primary photons for reconstruction as compared to a hardware-based approach.

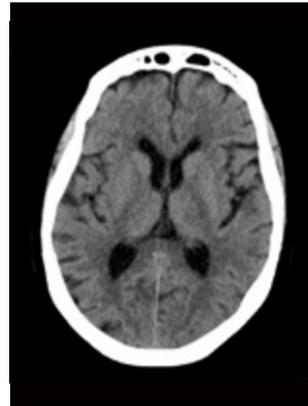

PUREVISION Optics – Delivering excellent brain image quality.



Eliminating the workflow challenges of MBIR* Integrated, easy to use, and fast



GENESIS Edition provides sharper image detail and lower patient dose with the world's first fully integrated MBIR solution.

FIRST utilizes forward projection iterations to deliver high-quality images with up to 82% dose reduction. A full volumetric reconstruction for routine clinical use can be obtained in just 3 minutes.

Following Toshiba's longstanding philosophy of minimizing dose while maintaining efficient clinical workflow, FIRST integrates seamlessly into your daily clinical practice.

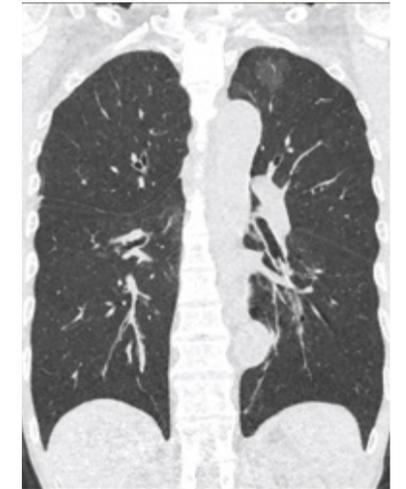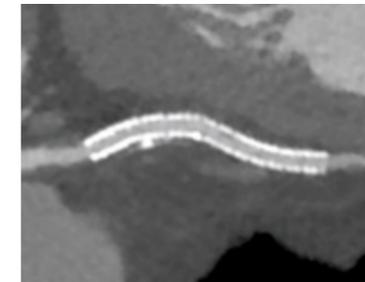

Chest scan performed with 0.2 mSv.

FIRST – Forward projected model-based Iterative Reconstruction SoluTion.
* Model-Based Iterative Reconstruction

The right application for a confident diagnosis

Automated, reliable, and robust


GENESIS Edition offers a comprehensive suite of Adaptive Diagnostic solutions to make complex exams easier and to improve diagnostic precision and reproducibility.

Subtraction CTA*




Superior visualization in CTA with true subtraction of bone and calcification.

Iodine Mapping*

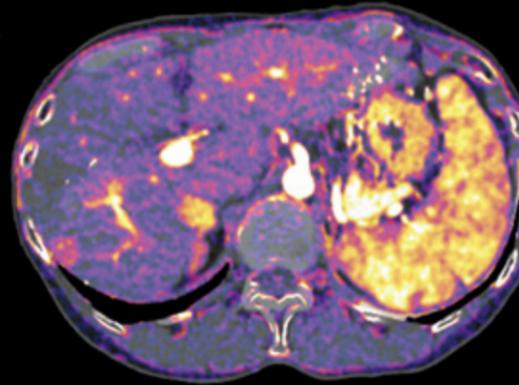

Clearly defined perfusion with color blood flow maps as a result of advanced registration and subtraction.

SURE Cardio™*

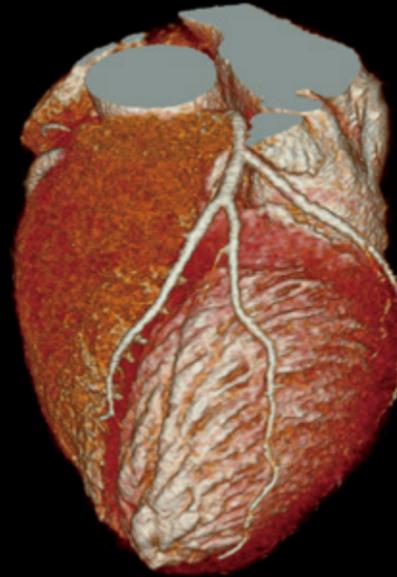


The robust solution for coronary imaging with ONE shot volume imaging and arrhythmia scanning.

TAVR*



Easily combined gated and non-gated acquisition for fast and low-dose TAVR exams.

SEMAR™

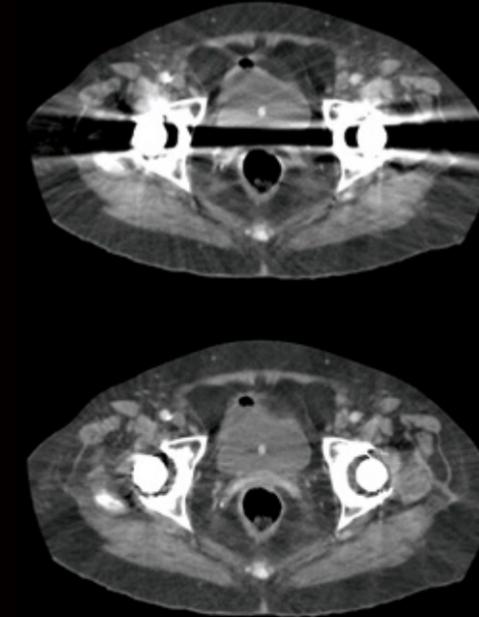

Improved visualization of bone and soft tissue-Single energy raw data based metal artifact reduction.

DE Tissue Characterization*





Tissue Characterization with easy-to-use Dual Energy scanning.

Dynamic Volume CT – Simply efficient

More than a decade of clinical partnerships with leading institutions sets Toshiba apart as the industry leader in dynamic volume CT. Together, we have developed new procedures for better patient care, automated workflows, and refined reconstruction technology to make the remarkable routine.


Aquilion ONE / GENESIS Edition – A new standard in CT based on 10 years of clinical know-how
There is no substitute for experience.

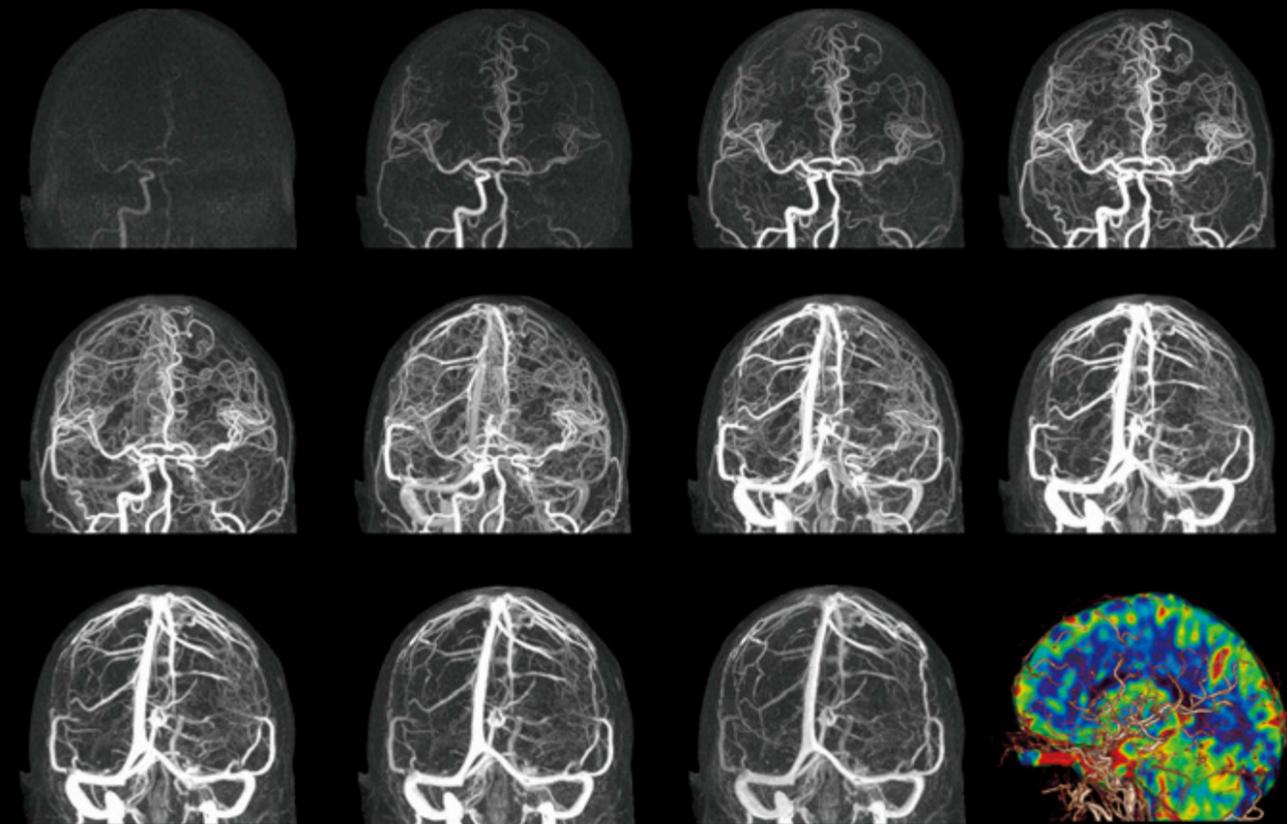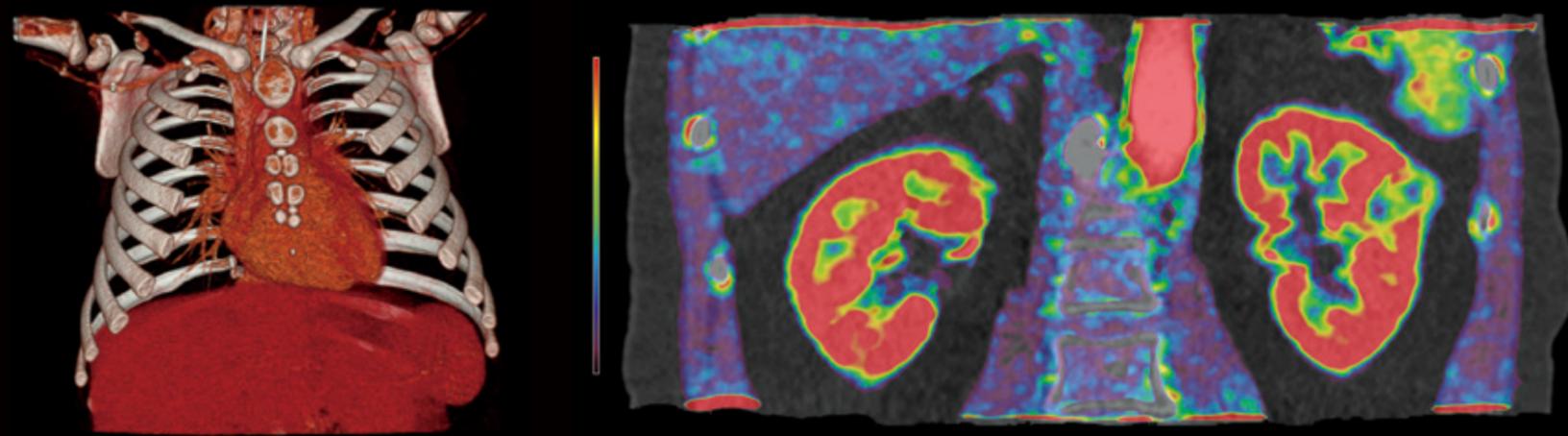




Patient-centric design

GENESIS Edition has been designed with a unique flared gantry, providing a calming, wide-open space for a better patient experience. The short bore is safer, with improved access to the patient from the front and rear of the gantry. During trauma and interventional procedures, patients can easily be cared for from the front and rear of the gantry.

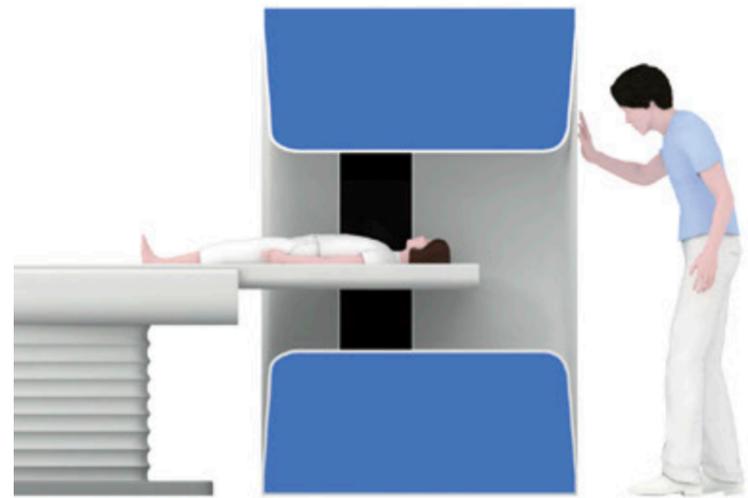

Claustrophobic tunnel design

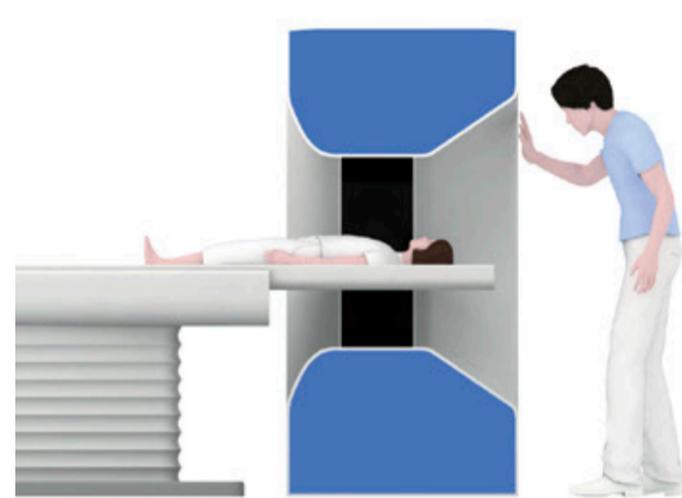

GENESIS Edition – Patient-friendly open flared design




The open flared design ensures that all patients, from the youngest to the largest, will remain at ease during scanning and provides unobstructed access at all times.

A couch designed for patient and technologist safety

Tech Assist Lateral Slide* ensures safety and comfort by providing a tool to mechanically move the patient to the correct position with the push of a button.

Tech Assist Lateral Slide reduces the risk of injury to the patient and the technologist.

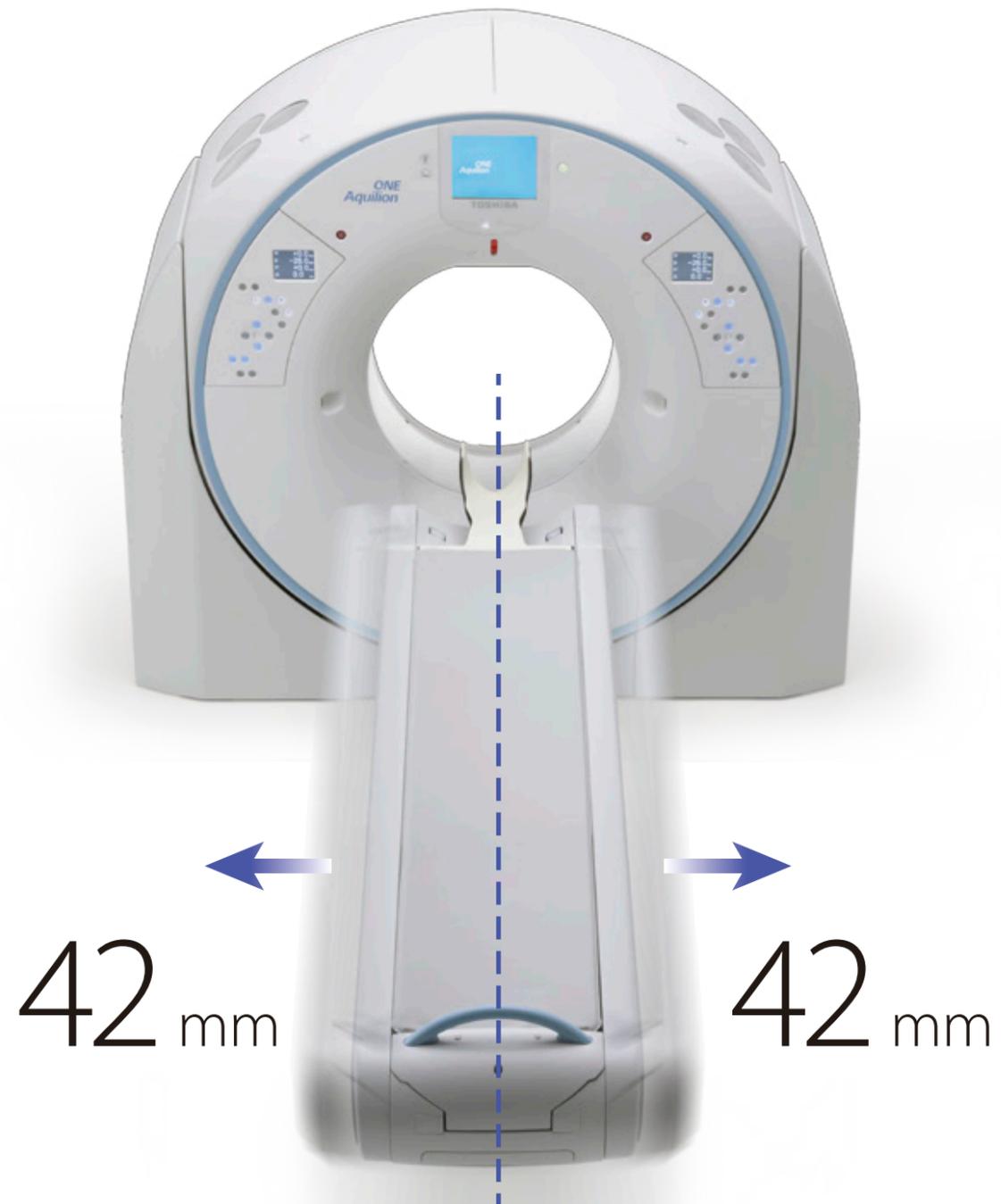

*Option

No need to accept a premium CT system without gantry tilt

Gantry tilt is a fundamental feature of CT systems, permitting angled scanning at your desired reading plane and avoiding direct exposure to radiosensitive organs.

Precision engineering equips GENESIS Edition with bidirectional gantry tilt. Highly advanced reconstruction technology overcomes the mathematical complexity of angled scanning for helical and volumetric acquisition, with no compromise in image quality.






Introducing the ease of X-ray in CT


Thanks to GENESIS Edition, CT scans can now be performed with the simplicity of conventional X-ray. Laser collimation* allows the field of view and scan range to be set directly on the gantry. Patients are positioned more comfortably. Examinations are performed faster and with reduced radiation dose.

*Option

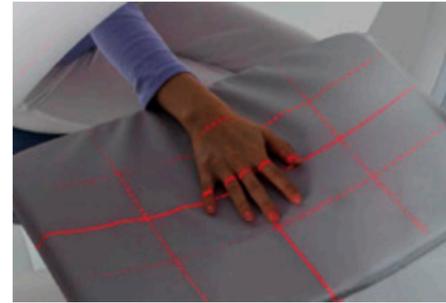



Transforming your workspace

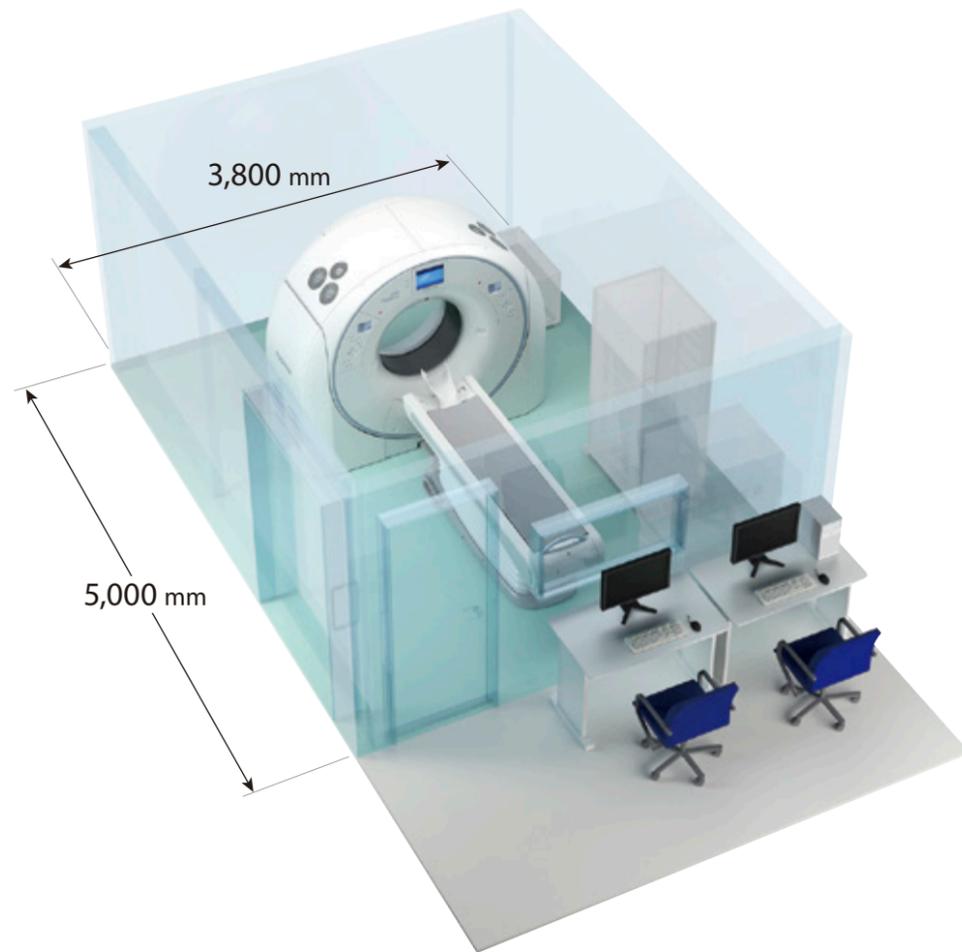

Economize on space, not on performance

GENESIS Edition is smaller, lighter, and requires less power than any other premium CT system. Designed for an installation space of just 19 m², GENESIS Edition can be installed in most existing CT rooms, avoiding costly renovations.

The compact design also provides more in-room space for trauma or interventional procedures.



19 m²
Installation Space



125 kVA
Power Capacity



2,360 kg
Gantry Weight





Making your workflow


GENESIS Edition helps to make complex exams easier, reducing dose and improving diagnostic accuracy and reproducibility. All steps from exam planning to reconstruction and postprocessing can be combined in the same protocol. So simply selecting from the wide range of pre- or user-defined protocols is all that's needed to achieve rapid and robust results.

Powerful Protocols


Patient Adaptive Exposure

Advanced Applications

AUTOMATED ▶▶▶

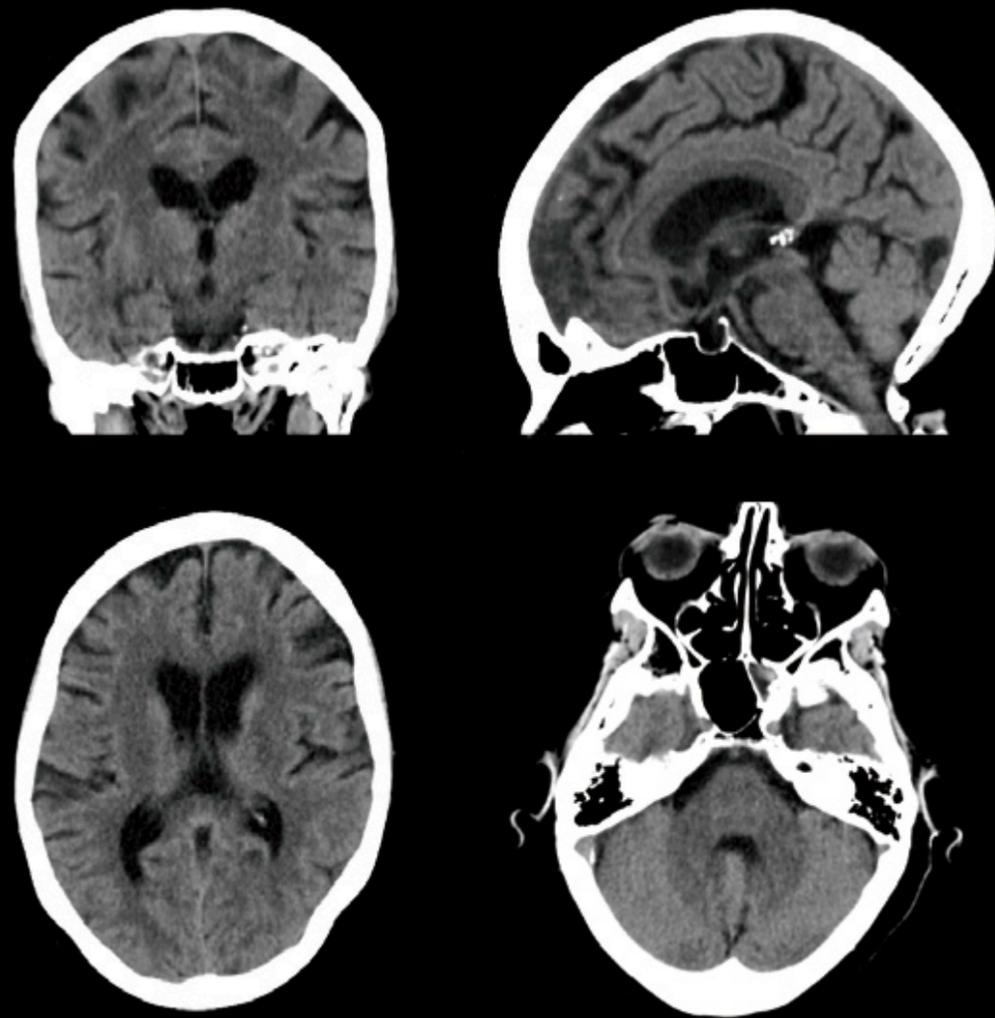





Superb brain image quality with clear grey-white matter differentiation and significantly reduced artifacts thanks to PUREVISION optics

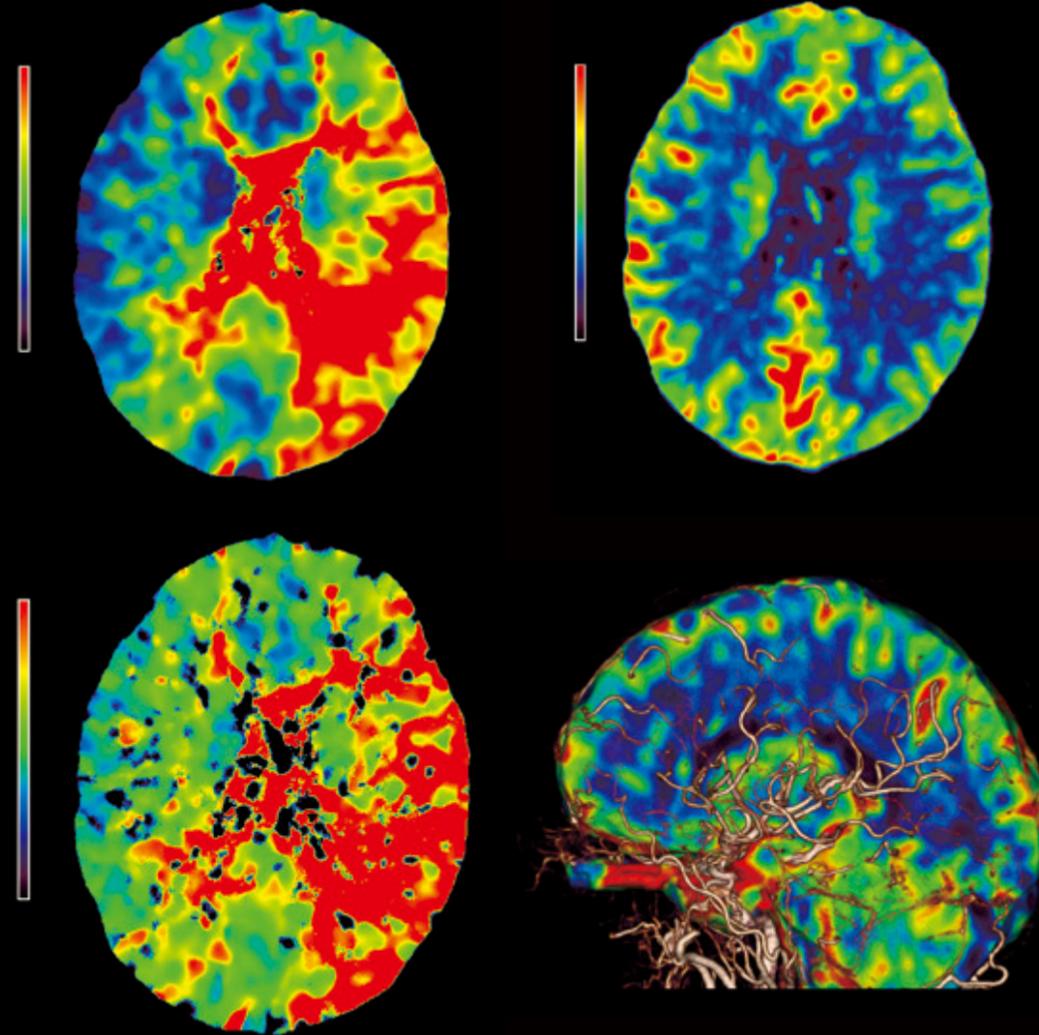

Whole-brain perfusion* with Dynamic Volume CT

*Option

Chest Imaging

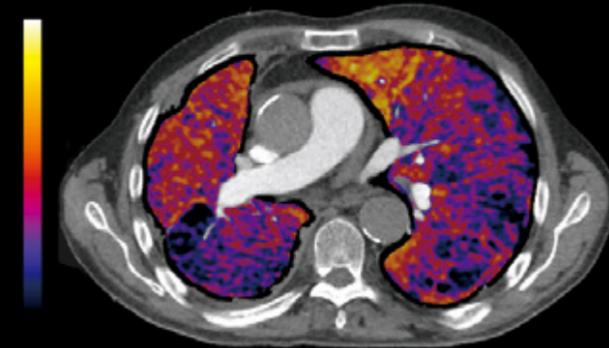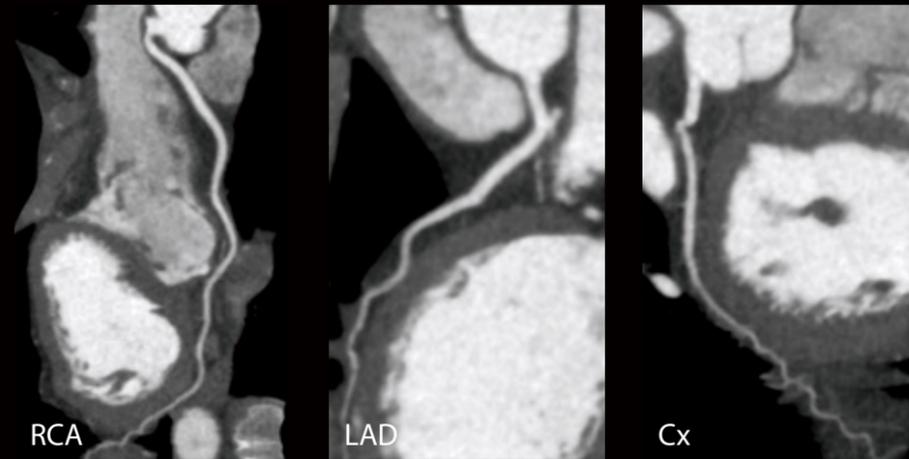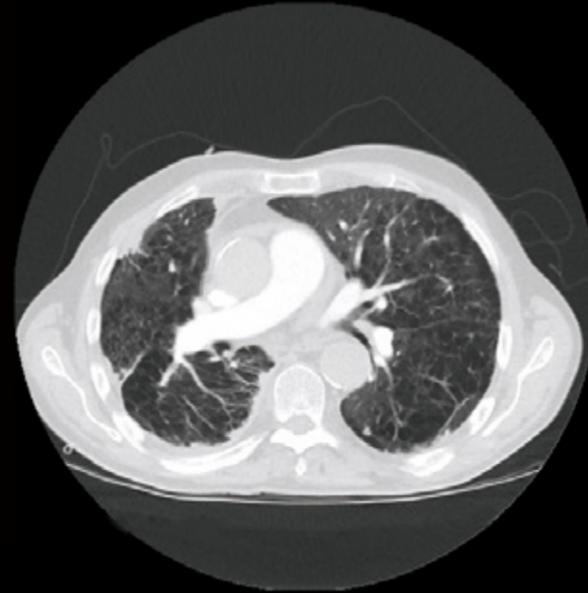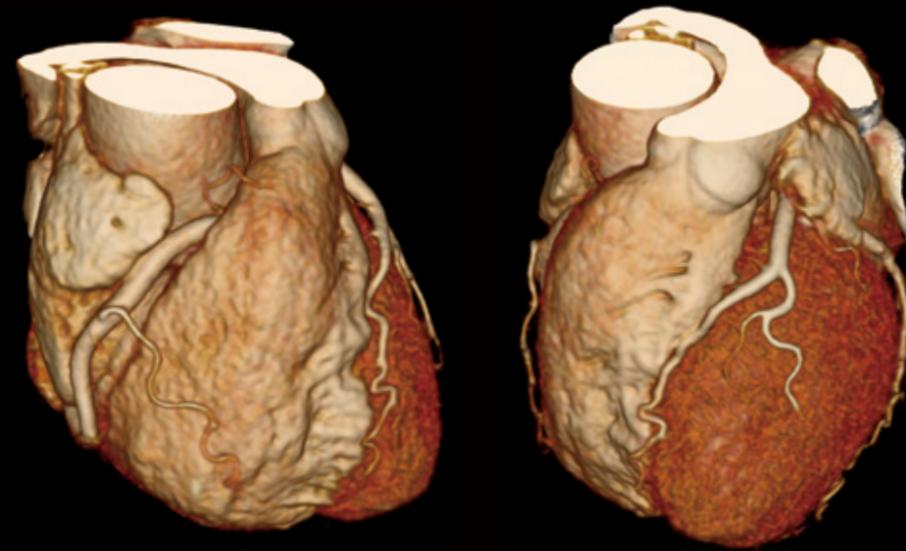


ONE beat cardiac CT* with high spatial resolution at only 0.26 mSv

Enhanced clinical confidence with ^{SURE}Subtraction™ Lung*

ONE volume chest scan of newborns: ultra-fast scan, ultra-low dose

*Option

Body Imaging

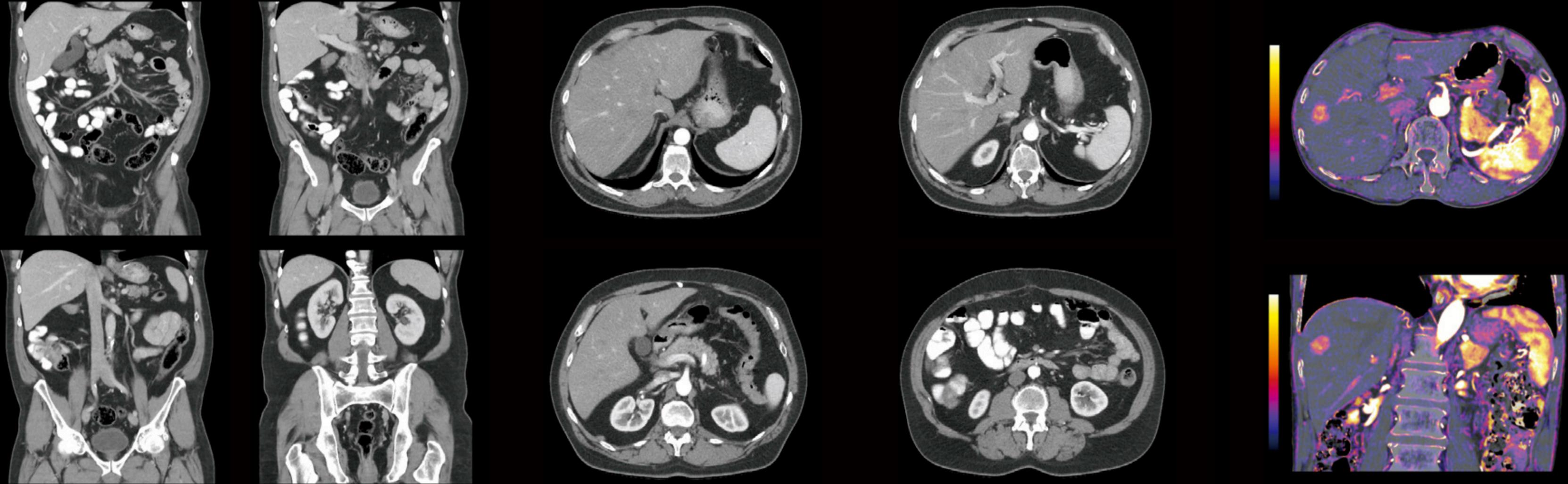

Excellent image detail and low contrast resolution in the abdomen

Increased lesion detectability with ^{SURE}Subtraction Iodine Mapping*




*Option



Original





SEMAR

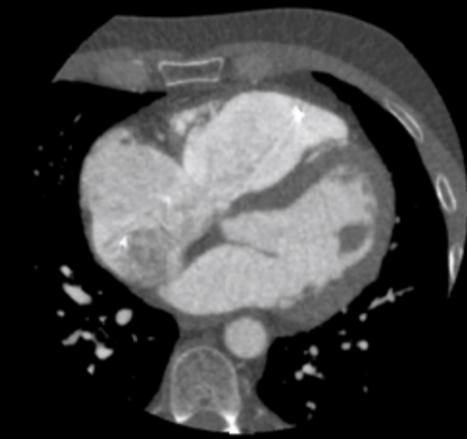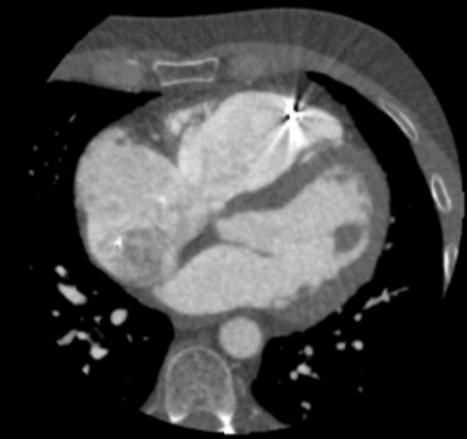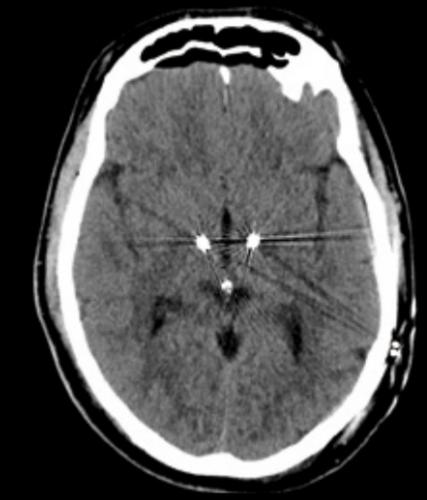

Original



SEMAR

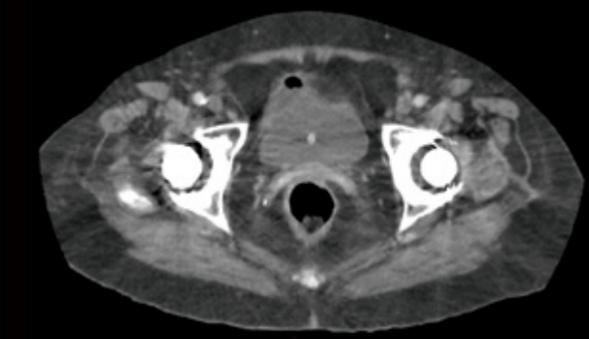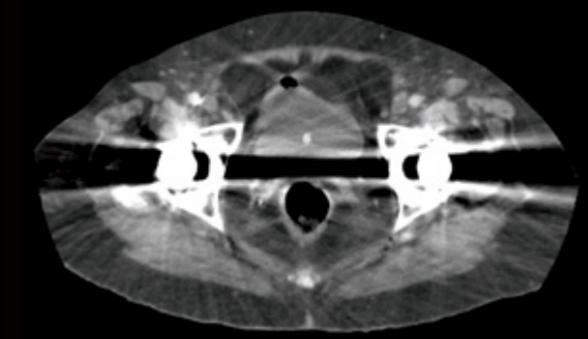



Main specifications		
Detector		PUREVISION detector
		320 rows, 0.5 mm
Gantry	Rotation time	Min. 0.275 s
	Bore size	78 cm
	Bore depth	38.7 cm
	Tilt	± 30°
Patient couch	Load	300 kg
	Max. scan range	150–200 cm
Reconstruction speed	Volume	5 s
	Helical	Max. 80 fps
Reconstruction	Iterative reconstruction	AIDR 3D* Enhanced
	MBIR	FIRST
Installation	Power capacity	125 kVA
	Space	Min. 19 m ² (short couch)

*Adaptive Iterative Dose Reduction 3D



GENESIS Edition Transforming CT

WARNING: Any reference to x-ray exposure, intravenous contrast dosage, and other medication is intended as a reference guideline only. The guidelines in this document do not substitute for the judgment of a healthcare provider. Each scan requires medical judgment by the healthcare provider about exposing the patient to ionizing radiation. Use the As Low As Reasonably Achievable (ALARA) radiation dose principle to balance factors such as the patient's condition, size and age; region to be imaged; and diagnostic task.

Disclaimer: In clinical practice, the use of the AIDR 3D and FIRST feature may reduce CT patient dose depending on the clinical task, patient size, anatomical location and clinical practice. A consultation with a radiologist and a physicist should be made to determine the appropriate dose to obtain diagnostic image quality for the particular clinical task. Due to local regulatory processes, this product may not be available in each country. Please contact your local Toshiba sales representative for the most current information.

TOSHIBA MEDICAL SYSTEMS CORPORATION

<http://www.toshibamedicalsystems.com>

©Toshiba Medical Systems Corporation 2016. All rights reserved.
Design and specifications subject to change without notice.
Model number: TSX-305A MCACT0289EAA 2016-07 TMSC/D

Toshiba Medical Systems Corporation meets internationally recognized standards for Quality Management System ISO 9001, ISO 13485.

Toshiba Medical Systems Corporation Nasu Operations meets the Environmental Management System standard ISO 14001.

Aquilion ONE, Aquilion ONE GENESIS, ^{SURE}Cardio, SEMAR and ^{SURE}Subtraction are trademarks of Toshiba Medical Systems Corporation.